



Designed to Last



Design

The VIDAS® range is greatly
Even with years of experie
So, after listening to your c
we've taken this already in



Unique values with a stylish-design puts VIDAS

Improved Design for Increased User Comfort

From your requirements we have advanced the system's design to further :

- **Streamline your workflow** by optimising test loading
- **Increase user-friendliness** for data entry and reporting patient results.

Timeless VIDAS Values

- Single-dose, ready-to-use reagents (1 patient = 1 test)
- All inclusive reagent kits
- System robustness (MTBF* > 452 days)
- Calibration every 14 days
- Automated barcode identification
- Just Load and Go

Flexibility, where and when you need it :

- Five independent 6-test sections
- Processing of different parameters simultaneously
- Large routine menu (over 80 parameters in Immunochemistry and Infectious Diseases)
- Packaging suited to different test volumes (30 or 60 tests)

Quality

- No cross-contamination
- Reference tests : HIV DUO Ultra, D-Dimer Exclusion™, Toxo Avidity...
- CVs highlighted as being excellent by numerous external quality controls

* Mean Time Between Failure

VIDAS Size:


58 cm (23 inches) (height)
105 cm (41.5 inches) (width)
68 cm (21.75 inches) (depth)

VIDAS weight:

65 kg (143 lbs)

Electrical requirements:

Voltage: 100-240 VAC
Electrical Consumption : 3-1.2 A
Frequency: 50/60Hz
Power: 280 VA maximum
Heat emission: approx.921 BTU/hr

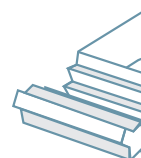


REAGENT

...d to Last

...y appreciated worldwide for its simplicity, flexibility and accuracy.
...nce, our continued success relies on you.
...urrent needs to further optimise workflow
...novative design one step further.

AS in a class of its own.

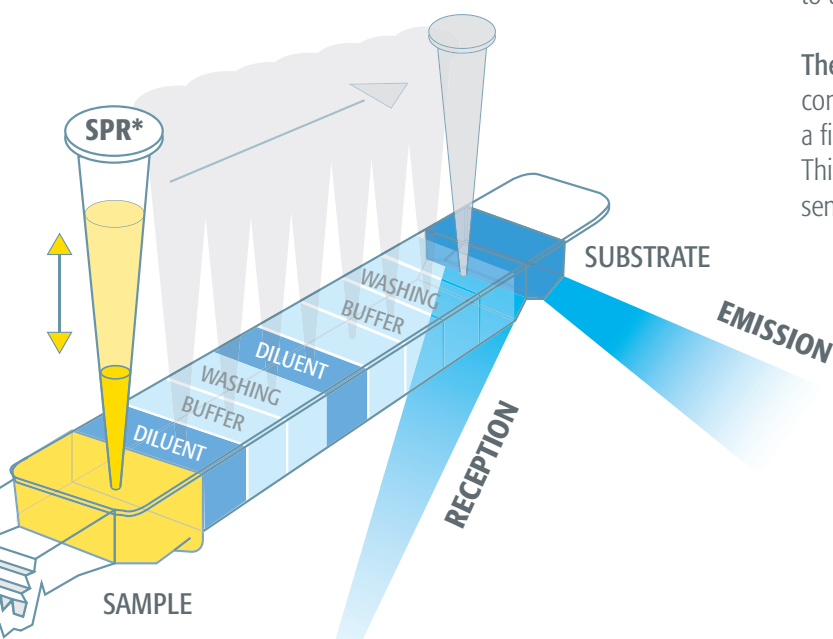

Simplicity makes the difference

The VIDAS test system uses a single-dose, ready-to-use reagent.

- 1• The **SPR**[®] acts as a Solid Phase Receptacle for the reaction and is coated with antigens or antibodies.
- 2• The **Strip** contains all the reagents required for the reaction.

The SPR acts as a pipetting device.

At each stage of the reaction, it aspirates the reagents in and out. This original concept prevents any inter-reagent or inter-sample contamination. In addition, the absence of tubing, syringes and needles reduces system maintenance to a minimum.

The VIDAS system uses the ELFA assay principle, combining the ELISA test method with a final blue fluorescent reading. This technology ensures excellent result sensitivity and specificity.

* SPR: Solid Phase Receptacle

VIDAS, Designed to last

Ease-of-use: Load and Go

Robustness: MTBF > 700 days

High Quality tests for reliable diagnosis

Cost effective solution

- **Single dose**
- **All inclusive kits**
- **Convenient packaging size**

VIDAS and mini VIDAS : adapted offer to your organisation.

VIDAS, a Customized Solution for your Daily Workload

Routine Testing

- Up to 10 parameters in one run
- Comprehensive panels (Thyroid, fertility, cardiac markers, tumor markers, HIV, hepatitis, ToRC...)

Complementary and Specialized Testing

- Confirmatory testing (HIV DUO, Toxo Avidity, DD Exclusion, Prolactin...)
- Specialized testing: Lyme, Measles, Mumps, H. pylori, CD A/B, Chlamydia, Allergy...)
- Complementary Testing: large menu (> 90 parameters)

Emergency Testing

- 24/7 available
- Rapid results in compliance with Healthcare guidelines
- Emergency panel: TNIU, CK-MB, Myoglobin, NT-proBNP, DD Exclusion, PCT, Digoxin, hCG, p24 Ag...



Z.A.C du BOIS MOUSSAY , 10 Avenue
du colonel ROL-TANGUY - 93240

STAINS FRANCE

Tel: +33149988672, +33953341938

Fax: +33149988672 - Mob:

+33678621481

socimed@socimed.com