

Floor Stand Series

Floor Stand Series - The Perfect Fit

The **Floor Stand Series** are large sized autoclaves built for **Hospitals and Medical Centers**. The 3 products in this series, with their flexible options and accessories provide plenty of diversity to allow you to mix and match to find the perfect fit for your organization.

The Floor Stand Series includes the **Pointer** (Fixed Pressure), the **SAP** (Pressure Control) and the **SAT** (Temp Control). The Pointer is the “Basic” version, followed by the “Mid Range” SAP and the “Full Function” SAT. The Pointer itself has two versions; the “Manual” version, which requires manual intervention during a session; and the “Auto” Pointer version, which is fully automated once a session is started.

- PLC Control
- Pre Vacuum
- Standby Jacket Warmer (SAT)
- Post Vacuum (standard on SAT, optional on SAP, Pointer)

For Hospitals and Med Centers


Everything you need for a right match

Product diversity, abundant options, complete accessories

FLOOR Stand


Mid-Range

SAP

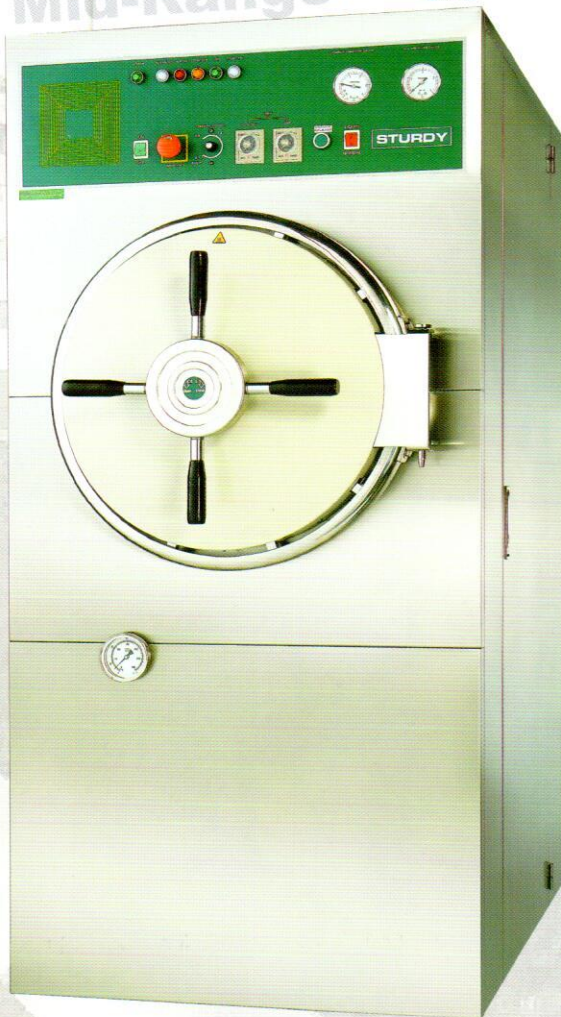

Full Function

SAT

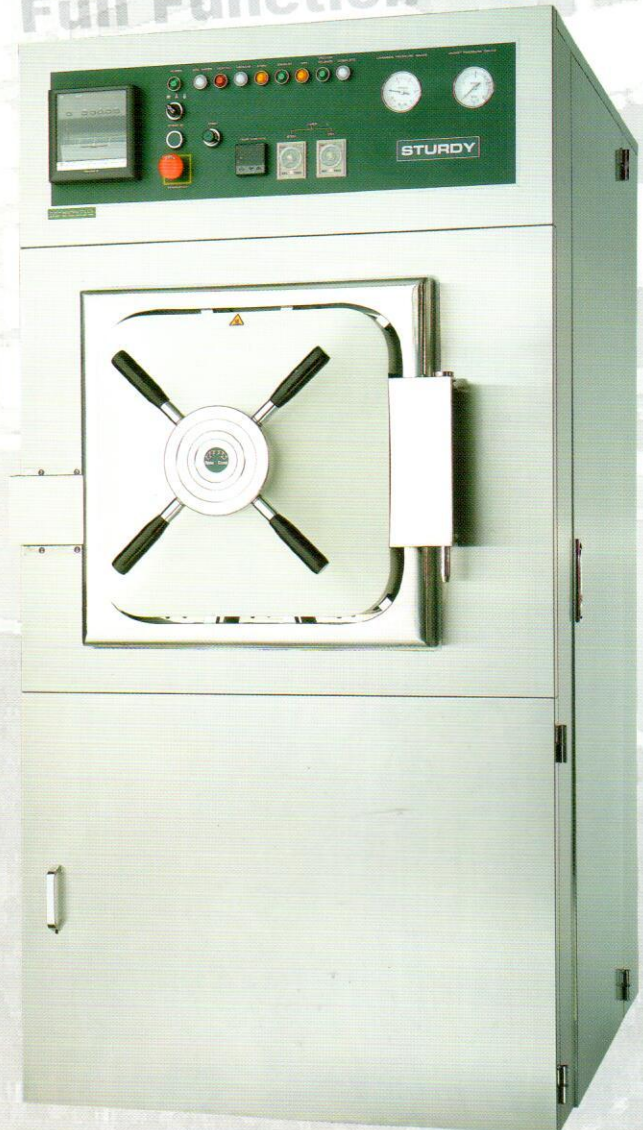

CE 0575 PED

For Hospitals and Med Centers

Floor Stand Series Feature Comparison

Model	SAT	SAP	Pointer-Auto	Pointer-Manual
Specification				
Control System	Digital Temp. PLC Control	Mechanical Pressure Control Optional (T) Digital Temp. Control	Mechanical Pressure Control	Mechanical Pressure Control
Program Selection Wrapped/Unwrapped/Liquid	Yes	Yes	N/A	N/A
Air Removal	Vacuum Pump	Gravity	Gravity	Gravity
Pre-Vacuum Function	Yes	N/A	N/A	N/A
Post(Dry)-Vacuum Function	Yes	Optional (V)	Optional (V)	N/A
Dry Function	Yes (Auto)	Yes (Auto)	Yes (Auto)	Yes (Manual)
1-point Temp. Recorder	Yes	Optional (X)	N/A	N/A
6-point Temp. Recorder	Optional (Y)	Optional (Y)	N/A	N/A
Double Door	Optional (D)	Optional (D)	N/A	N/A
Elec. Vertically sliding door (up-down)	Optional (F) - Square Chamber Only	Optional (F) - Square Chamber Only	N/A	N/A
Fill water system* (water pressure of 2 bars required)	Auto	Auto	Manual Optional (B) for Auto	Manual
Chamber/Jacket Pressure display	Gauge	Gauge	Gauge	Gauge
Stand-by Function (Switch)	Yes	N/A	N/A	N/A
Emergency/Reset Switch	Yes	Yes	N/A	N/A
Add-Water Pump	Optional (W)	Optional (W)	Optional (BW)	N/A
Water Level Protection	Yes	Yes	Yes	Optional (L)
Pressure Door Auto-Lock	Yes	Optional (K)	Optional (K)	Optional (K)
Life Time	7 Year	7 Year	7 Year	7 Year

Option Codes (e.g. Optional (V)) listed in the table above are only reference codes used by other documentation.


B: AUTO-ADD WATER PLUMBING (water pressure of 2 bars required)
D: DOUBLE DOOR DEVICE
F: ELEC. VERTICALLY SLIDING DOOR (SQUARE CHAMBER ONLY)
K: PRESSURE DOOR AUTO-LOCK DEVICE
L: WATER LEVEL CONTROLLER (LOW WATER INDICATOR)

T: DIGITAL DISPLAYED TEMP. CONTROLLER
V: DRY VACUUM PUMP
W: ADD WATER PUMP
X: 1-POINT TEMP. RECORDER
Y: 6-POINT TEMP. RECORDER

Floor Stand Series (continued)

Floor Stand Series Accessories

Square Chamber (SAT, SAP)

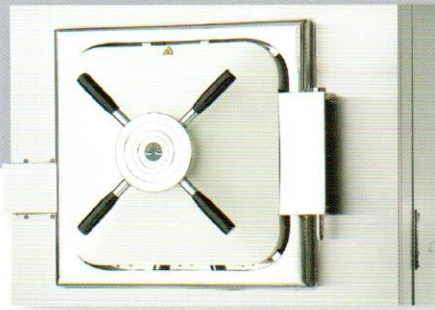

Basket




Tray



Plate Set



Round Chamber (SAT, SAP, Pointer)

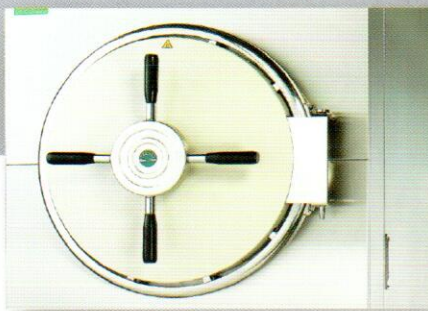

Basket



Tray



Optional
Trolley
(SAT, SAP
Pointer)



Temperature Recorder for SAT (standard) and SAP (optional)



Optional 6-pt Sensor

The Temp Recorder comes standard with a 1-point Temp Sensor. 6-point Temp sensors can optionally be installed in both the SAT and SAP.